DESIGNING GOOD LIFE

Shiatsu Massage TOTAL CONTROL

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Safety instructions

1. Read the instruction manual carefully before use, and save it for future use.
2. Plug the appliance to a household outlet of 220-240V only, and use the appliance only for the intended use.
3. Always unplug the appliance from the mains supply when disorder is experienced, when not in use, when parts are applied or removed from the appliance and before cleaning and maintenance.
4. The appliance must not be immersed into water or other liquids.
5. Persons with blood circulation problems or reduced sensibility, e.g. diabetics, should consult a physician before using the massage cushion. Persons with pacemakers should not use the massage cushion. Pregnant women should consult a physician before using the massage cushion.
6. This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they have been given supervision or instruction concerning the use of the appliance in a safe way by a person responsible for their safety and they understand the hazards involved.
7. Children should be supervised to ensure that they do not play with the appliance. Children are not always able to understand potential risks. Teach children a responsible handling of electrical appliances.
8. The appliance becomes hot during use. People who are insensitive to heat should use the device with caution.

9. Cleaning and user maintenance shall not be made by children, unless they are over the age of 8 and supervised.
10. The appliance should only be used in an upright position.
11. Do not use the appliance in bed.
12. Do not use the appliance when sleeping.
13. The appliance should never be left unattended when plugged in or when operating.
14. Do not keep the appliance next to a naked flame.
15. Do not smoke when using the massage cushion.
16. Always examine appliance, cord and plug for damage before use. If the supply cord is damaged, it must be replaced by the manufacturer, or his service agent or a similarly qualified person in order to avoid a hazard.
17. For additional protection the installation of a residual current device (RCD) with a rated residual operating current not exceeding 30 mA is recommended. Ask your installer for advice.
18. The appliance is for domestic use only, and must not be used professionally or medically.
19. If the appliance is used for other purposes than intended, or is it not handled in accordance with the instruction manual, the full responsibility for any consequences will rest with the user. Any damages to the product or other things are not covered by the warranty.

Before use

IMPORTANT! The massager comes with a screw located at the back of the unit to protect the massager mechanism due to shipment, which must be removed before the first use to allow the shiatsu feature to function. Use the provided Allen key.

Then properly dispose of the screw. **WARNING!** Failure to remove screw may cause permanent damage to the massage cushion.

Error: all remote control LEDs flash for a few seconds.

Cause: massager overloaded.

Solution: ensure the allen key bolt/screw has been removed from back of the massage cushion (using the allen key included). Should the condition continue after this, please contact the retailer.

Connect the plug to a wall socket. Choose between the following programmes from the control panel:

1. DEMO

A brief demonstration of all programmes.

2. Adjust the position of the massage heads

3. SHIATSU

Shiatsu massage for the whole back.

4. SHIATSU

Shiatsu massage for the upper back.

5. SHIATSU

Shiatsu massage for the lower back.

6. SHIATSU spot massage

7. Vibrating massage for the thighs.

8. Switch ON/OFF

9. ROLLING

Rolling massage for the whole back.

10. ROLLING

Rolling massage for the upper back.

11. ROLLING

Rolling massage for lower back.

12. Button for adjusting the width between the massage heads.

13. Heat function

Important! Persons with blood circulation problems or reduced sensibility, e.g. diabetics or persons with neck injury, should consult a physician before using the massage cushion. Persons with pacemakers should not use the massage cushion. Pregnant women should consult a physician before using the massage cushion. The body massage does not compensate medical treatment. At continuous pain and ache you should always contact medical expertise.

For use at home or office

OBH Nordica's Shiatsu massage cushion incorporates a unique strapping system which allows you to fasten it to almost any chair by tightening the Velcro fastening around the backrest.

As well as the strapping system to secure the cushion to the back of most chairs, the massager also has an adjustable seat which enables use on most soft arm chairs and sofas. Simply adjust the strap length using the buckle under the seat to achieve the desired angle to fit your chair/sofa and the cushion can be used effectively at virtually any angle. By all means, you can also just put the cushion in a chair or sofa if you prefer. The massage cushion is not intended for use in cars and may not be used lying down.

Please note: The massager is not recommended for use on wooden surfaces as the slipper may damage the wood. Caution is also recommended when using on upholstered furniture. Use some kind of thicker textile as protection.

Switch (8)

Connect the appliance. To turn on the massage cushion functions, press the switch. The pilot lamp lights. To switch off the massager, simply press the button again. The pilot lamp will blink while the dual massage mechanism returns to its lowest point and then it will turn off.

DEMO (1)

For a brief demonstration of the features of the Shiatsu massage cushion, press the DEMO button. The unit will complete one cycle combining Shiatsu, rolling and vibration massage. Once complete, the unit automatically turns off. The LED light will blink to indicate the active functions.

Adjusting the massage heads (2)

To customize the massage to your height, press the button at the highest desired position. The massage heads will not travel higher than that position for the duration of the massage. To deselect the feature press the button a second time. The massage will return to its full travel.

SHIATSU massage (3, 4, 5)

Shiatsu is a Japanese massage technique which means a deep kneading circular massage. Choose from 3 preset programs where the Shiatsu massage travels up and down in a specific area of the back. To select a program, simply press the button and the pilot lamp will illuminate. To deselect, press the button again and choose another program.

SHIATSU spot massage (6)

While the unit is in Shiatsu massage mode, press either up or down button to activate the spot massage. You may then adjust the mechanism position by pressing corresponding button until the desired position is reached. When you release the button the massager stops at desired setting.

Vibrating massage (7)

Adds vibrating massage to upper legs/thighs.

ROLLING message (9, 10, 11)

Rolling is like two hands or thumbs working up and down along your spine. Choose from 3 preset programs where the Rolling massage travels up and down in a specific area of the back. Whole back, upper back and lower back. To select a program, simply press the button and the pilot lamp will illuminate. To deselect, press the button again or choose another program.

Adjustment of width between the massage heads (12)

While the unit is in ROLLING massage mode, one can adjust the distance between the massage heads. This adjustment changes the area that is massaged along your spine. With each press, the massage heads will move to the next position.

Heat function (13)

Activate the heat function by pressing the button and the corresponding light will illuminate and the massage heads will glow red. To deselect press the button again. Note that the heat function only can be activated after choosing a massage program.

CAUTION!

To avoid pinching do not lean on the dual massage mechanism in the cushion when adjusting your body position. Do not jam or force fingers or hands in the moving massage mechanism.

The massage cushion Shiatsu and Rolling massage is only to be used at the back/neck.

If the fabric or material covering the massage mechanism becomes torn or damaged, discontinue use.

Only gentle force should be exerted against the unit in order to eliminate any/all risk of injury and to avoid overloading the motor.


Soften the massage force by placing a towel between yourself and the unit.

Important! The cushion will auto shut off after 15 minutes. This because the massage cushion should not be used for more than 15 minutes at a time. After a 15 minutes treatment we recommend a rest for at least 15 minutes Note! The moving dual massage mechanism always "parks" or finishes in its lowest position.

Cleaning and maintenance

Store the massager in its carry bag (see instructions below) or in the box in a safe dry place.

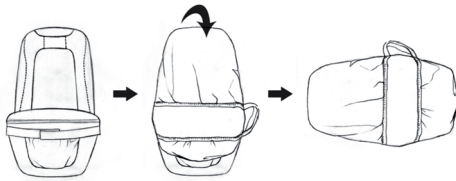Avoid contact with sharp edges or pointed object that might cut or puncture the fabric surface. To avoid breakage do not wrap the power cord around the unit. Pull the plug before cleaning. The massager can be cleaned with a damp cloth. After dry with a towel. Water or other liquids can not come in contact with the massager. Never use any kind of detergent when cleaning.

The removable back flap can be washed. Simply follow the washing instructions:

Disposal of the appliance

Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.

Warranty

Claims according to current law. Cash receipt with purchase date shall be enclosed with any claims. When raising a claim the appliance must be returned to the store of purchase.

OBH Nordica AB

Löfströms Allé 5
172 66 Sundbyberg
Tel: +46 08-629 25 00
www.obhnordica.se

OBH Nordica 6092
220-240 V AC, 50Hz
35 watt
Uptime: 15 min per treatment.

With reservation for running changes.

SE/ES/6092/UVN/020915